


comationTM

[changing Lives With Technology]

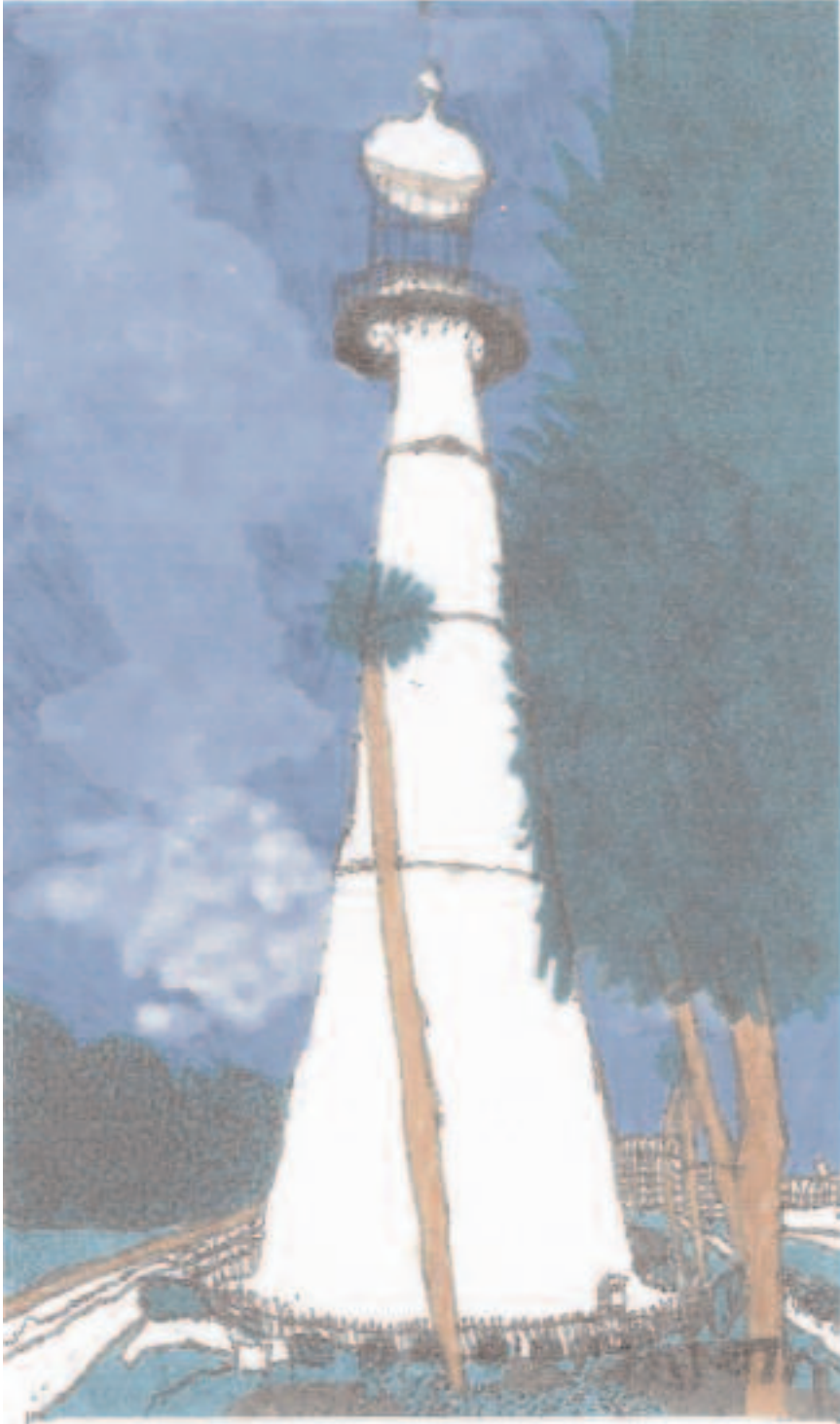

*Biloxi Lighthouse by Salome McCaffrey

Meet Salome...

Salome Alexandria McCaffrey

Born 12-25-1983 Early Christmas morning.

At 6 months she was ahead of all the other children her age. At 9 months she was back to a newborn stage. Doctors said she would not make it a year. At 1 year old she received plasmafloresis. At 18 months she was taped to be on the Jerry Lewis MDA television marathon, although she never got to see it. Hurricane Elena came through town and took hers lungs with it. She was ventelated and trached. She is 20 and still going strong. Thank God.

In the fall of 2001 The Make a Wish Foundation arranged for Salome to have her wish which was to meet Rucker Ashmore, Lisa Rotelli, and the company that helped Salome to be independent & free. She thought he ought to know how important their work was for her and others like her.



Contents

2006 CATALOG



Communication 3
Rolltalk communication aids, environmental control, software, and accessories



Mechanical Switches 13
All shapes sizes and activations.



Sensors 9
No touch sensors and fiber optic switches.



Mounting Hardware 15
Mounting hardware for all switches.



Power Options 10
Power options for all sensors.



Computer Access 16
Mouse emulators hardwired and wireless.



Packages 12
All in one packages and complete systems.



Fun N Games 17
Toys and Games for kids to use with switches.

THE POWER TO EFFECTIVELY COMMUNICATE....

has given Bobby the ability to control a few things in his life. Bobby has been able to let us know when he is getting sick rather than having to wait for his symptoms to become visually apparent to us. He has also enjoyed taking control of the television and stereo not to mention playing with a remote control car and dog.



ROLLtalk



The Rolltalk Mini

The all encompassing software allows the client to move fluently through all layers of access. This system allows a client to completely and seamlessly maneuver through all activities such as writing and printing a document from a communication page, sending and receiving text messages and taking a picture or video and adding text from your communication pages. With an endless amount of IR memory, this system can have any amount of preset or direct selected TV channels, stereos and DVDs while still having communication on the same page. The possibilities are endless. Through the Rolltalk you can truly be independent in your environment, no more multiple layers and changing screens or formats to print, call, and communicate.

Rolltalk is an exceptional combination of hardware and software that is a totally integrated system.

- Mira Shah MSCCC-SLP Pediatrics Speech Pathologist in Portland, OR

Part# Rolltalk Mini \$8,500.00



The Rolltalk Lite

This system has the power of a mini and a nice compact size. The Rolltalk Lite's size makes it very transportable and the screen clarity is unmatched. It has unlimited IR functions, cell phone sms, and communication with the power of Rolltalk software.

Part# Rolltalk Lite \$6200.00



The Rolltalk Micro

This is an amazing hand held device with quite a voice. It is an augmentative communication device that allows you to communicate, and have unlimited environmental control. It comes with the amazing RollTalk software so creating custom pages is a breeze. It can be accessed through switch scanning or touch screen.

Part# Rolltalk Micro \$4000.00



The RollTalk Software PC Version

The all encompassing software allows a client to move fluently through all layers of access. With the Rolltalk software you can build custom pages for communication and multiple layers all on one screen. It has the power of the Rolltalk workshop without IR capabilities.

Part# Rolltalk PC Software \$2,500.00

ROLLtalk Maxi

The Maxi



Rolltalk Drive Module

The Maxi Rolltalk combines the Rolltalk drive module with the Rolltalk Maxi tablet PC and the Rolltalk Workshop software. This system allows you to drive a power wheelchair, access all seating functions, communicate, access a cell phone, sms text message, take photos, take videos, control your TV, Stereo, DVD, and access all other functions of the Rolltalk software through scanning and a single switch. With the addition of any head mouse, you can control the functions of the maxi Rolltalk with head movement and a single switch.

With the addition of an eye gaze system such as Quick Glance 3 from Eye Tech or Erica from Eye Response you can control the functions of the Maxi Rolltalk through eye movement ; with or without a switch.

(patent pending)

Rolltalk Drive Module

Part# Rolltalk Drive Module \$8,000.00

Rolltalk Tablet PC & Software

Part# Rolltalk Mini-DM \$8,500.00



Driving With Eye Gaze

Driving with eye gaze is an incredible new technology developed by Rolltalk.

This new technology will enhance the lives of many people. However, it is not as simple as portrayed on television. Eye gaze is best used indoors due to changing lighting factors. It is also best used with someone with very little head movement and it is highly recommended that you use a switch with the eye gaze for safety. Rolltalk is so functional that you can set the screen to work with eye gaze inside and scanning while outside or in bright light. The scanner requires a switch to be activated. You must specify wheelchair brand and type of electronics when ordering because not all wheelchairs are compatible. This new technology is a major advancement for independence especially for clients with ALS. Thanks to Rolltalk it is here!!
(available only through Com'ation.)

(patent pending)



ROLLtalk Access Options



An eye gaze system such as Quick Glance 3 from Eye Tech can be used in conjunction with a Rolltalk maxi. You can access all functions of the Rolltalk software including driving a power wheelchair and all of the seat functions of that chair.

Part# Quick Glance 3 \$6,980.00



An eye gaze system such as Erica Eye Response can be used in conjunction with a Rolltalk maxi. You can access all functions of the Rolltalk software including driving a power wheelchair and all of the seat functions of that chair.

Part# Erica Eye Response \$7,900.00



Touch Screen

This comes standard on all Rolltalk devices from the Micro to the Maxi.



Head Mouse

Any head mouse like this one available from Com ation, will plug directly into any of the four Rolltalk USB ports

Part# Head Mouse \$1,190.00



ME500

The mouse emulator is a great way to access a Rolltalk. The ME500 can be set up as a three or five switch mouse. This version of the mouse emulator must be directly connected to the Rolltalk through the USB port and requires no software.

Part# ME500 \$365.00



Other Methods Of Access

Single switch scanning, dual switch scanning, Morse Code are just some of the other access methods that are built in.

— — MORSE
 — . — — . . CODE

ROLLtalk Accessories



Web Cam

This web cam is an exciting accessory to Rolltalk. It allows the client to take photos and video and have a live video window on their screen. It is a great way to promote communication and explain activities or use for a wheelchair backup camera.
(Available through Com'ation)

Part# Rolltalk Web Cam \$150.00



Blue Tooth USB Adaptor

Blue Tooth cell phones like this one, come ready and enabled to send and receive text messages and cell phone calls directly through Rolltalk seamlessly.
(Blue Tooth USB adapter is available through Com'ation)

Part# Rolltalk USB-BT \$100.00



Blue Tooth Printer

This example of a commercially available Blue Tooth printer would allow the client with their Rolltalk to print remotely from anywhere within 50 feet of the printer.
(Not available through Com'ation)



Robo-Sapien

Any infrared toy or game like Robo-Sapien can be easily programmed to function through Rolltalk's trainable infrared system.
(Not available through Com'ation)



Tash IR Telephone

- * Meeting all your communication needs with our advanced LCD, infrared telephone
 - * Direct dial any number
 - * Store 60 names and numbers in directory
 - * Head set compatible
- (Available through Com'ation)


Part# Tash IR-Phone \$850.00

THE POWER of COMPUTER ACCESS....

I like the mouse emulator because it let s me do my school work on my computer by myself
- Caitlyn



LANGUAGE Generating Software

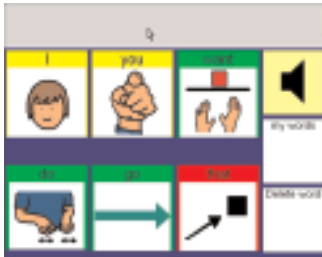

Talkabout is a new language program designed to keep the conversation going. With six locations, a variety of core words come to the user so that she or he may construct phrases and sentences in an error free manner. Choose one word and the locations change to predict six new logical word choices. Color-coded word categories make it easy to teach language while the user actively participates with generative and novel word combinations.



Coming soon generative language software built specifically with the Rolltalk workshop.



THE TRAIN CONDUCTORS...

Stormy and Ashley

Train conductors, Stormy and Ashley, are gearing up to catch the train as it comes around the track; timing is everything! Anticipating their photo, Stormy and Ashley wait and take the picture before the train gets away. All aboard...Whistles...Bells...Forward...Reverse... The train reacts to both conductors, Ashley and Stormy, as they activate various commands. The Rolltalk screen can be customized for multiple needs, direct selection or switch input!



Gatlin Flip Up Wheelchair Mount

This flip up mount is an alternative to a conventional communication mount. It mounts to the arm rest and pulls up and flips away to allow for transfers. If you add the tube clamp mount with a quick pull pin for quick release any device is tranportable.

Part# Gatlin Flip Up \$430.00



Mini Rolltalk Sleeve

Mini Rolltalk sleeve with heavy duty table mount easily transfers from wheelchair to the heavy duty table mount with the pull of a pin. The sleeve helps protect the Rolltalk from impact.

Part# Rolltalk Sleeve \$200.00



Part# RTC

Part# RTC-H

Tube Clamp Mount

In conjunction with the Rolltalk sleeve the tube clamp mount will allow a quick release to take the device from the chair to a table mount or over bed mount with the pull of a pin. If you add the turn handle you can easily adjust the angle of your device without tools.

Part# RTC \$135.00

Part# RTC-H \$155.00



Mini Rolltalk Wire Stand

This portable, lightweight table mount is great for on the go. You can fold it up for transport and place it on any flat surface to use your Rolltalk.

Part# Rolltalk Wire Stand \$85.00

Rolling BedSide Table Mount

This rolling height adjustable table is a perfect mount for a client to use their Rolltalk in bed. It is durably clamped on to the table so it is stable and will easily quick release to be transferred back to the powerchair.

Part# Rolltalk Bedside Mount \$450.00

With Tube Clamp Mount \$585.00



NO TOUCH NO PRESSURE *One Power Option Required

*ALL ELECTRONIC SWITCHES REQUIRE A POWER OPTION FROM Pg. 10 OR 11 TO FUNCTION.

Small Silver Round Proximity Sensor* 2" in Diameter x 3/8 Tall



This single function switch requires no touch or pressure. This proximity switch can detect something that has an electronic capacitance, such as a hand or finger, up to 1/4" away. It requires a power source from the power pages 10 or 11.

Part# PRX-101 \$165.00

*** MOISTURE WITHIN THE ACTIVATION RANGE WILL CAUSE SWITCH ACTIVATION. DO NOT COVER THIS SENSOR WITH ANY MATERIAL.

Small Silver Round Proximity Sensor* 1 3/16" in Diameter x 1/8 Tall



This sensor is like the PRX-101 but is smaller in diameter and thinner. This single function switch requires no touch or pressure. This proximity switch can detect something that has an electronic capacitance, such as a hand or finger, up to 1/4" away. It requires a power source from the power pages 10 or 11.

Part# PRX-102 \$165.00

*** MOISTURE WITHIN THE ACTIVATION RANGE WILL CAUSE SWITCH ACTIVATION. DO NOT COVER THIS SENSOR WITH ANY MATERIAL.

Non-Adjustable Proximity Sensor* 2" long- 3/8" thick- 3/4" wide



This single function switch requires no touch or pressure. This proximity switch can detect something that has an electronic capacitance, such as a hand or finger, up to 1/4" away. This sensor ignores non-capacitive materials and is ideal to mount in a headrest or tray. It requires a power source from the power pages 10 or 11.

Part# PRX-103 \$165.00

*** MOISTURE WITHIN THE ACTIVATION RANGE WILL CAUSE SWITCH ACTIVATION.

Adjustable Proximity Sensor* 1 3/8 wide 2 3/8 long 9/16 Tall



This single function switch requires no touch or pressure. This proximity switch can detect something that has an electronic capacitance, such as a hand or finger, up to 1/4" away. This sensor ignores non-capacitive materials and is ideal to mount in a headrest or tray. It requires a power source from the power pages 10 or 11.

Part# PRX-104 \$165.00

*** MOISTURE WITHIN THE ACTIVATION RANGE WILL CAUSE SWITCH ACTIVATION.

Large Fiber Optic Switch*



1/4 - Diameter x 5/8 - long

This fiber optic switch, with a round target area the size of a pencil, has an adjustment range of touch to 5 inches. This switch can be used in either the N.O. or N.C. mode. In the N.O. mode switch activation occurs when you interrupt or come into the beam of light. In the N.C. mode switch activation occurs when you come out of the beam of light. (NO-Normally Open / NC- Normally Closed)

Part# FOS-301 small fiber optic cable \$255.00

Part# FOS-302 large fiber optic cable \$255.00

*** DIRECT SUN LIGHT OR HIGH INTENSITY DIRECT LIGHTING MAY BLIND THIS SENSOR. BLINDING WILL NOT CAUSE FALSE ACTIVATION OF THE SWITCH; RATHER THE SWITCH WILL FAIL TO OPERATE UNTIL THE DIRECT LIGHTING IS REMOVED.



Power Options for Single Switch Sensors

*Electronic Switches Require Power from this Page



SWITCHES NOT INCLUDED / MUST PICK ONE FROM Pg. 9

Power 100

This portable power pack is a all in one box. It is rechargeable and has an option for a auditory cue that sounds whenever the switch is activated. It will last up to 72 hours of constant switch use and comes with a charger.

Part# PO-100 \$235.00




SWITCHES NOT INCLUDED / MUST PICK ONE FROM Pg. 9

Power 110+S

This portable power pack is ideal for long term sensor use. This battery will last up to 120 hours on single switch use before having to be charged. It comes with a wall transformer for recharging.

Part# PO-110+S \$200.00

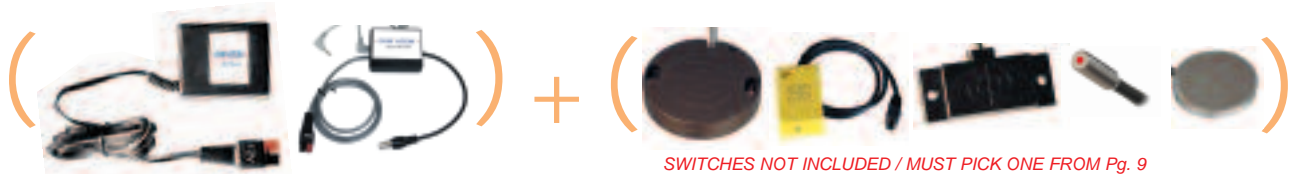

SWITCHES NOT INCLUDED / MUST PICK ONE FROM Pg. 9

Power 112+S Wall Plug

This wall plug is ideal for a stationary setting and will allow you to simply plug one of the electronic sensor outputs, (shown above) into the wall.

Part# PO-112+S \$75.00



SWITCHES NOT INCLUDED / MUST PICK ONE FROM Pg. 9

Power 130+S Powerchair Adapter

This is a portable system that will allow you to run a sensor from a power wheelchair.

Part# PO-130+S \$165.00

Chair Specific



Power Options For Multiple Switch Sensors

*Electronic Switches Require Power from this Page



SWITCHES NOT INCLUDED / MUST PICK ONE FROM Pg. 9

5 - Sensor Output/ Power 110

The five switch adapter will allow you to connect up to five sensors through one hub. This is a portable package and will allow you to go anywhere.

Part# SO-5PR+110 \$350.00



SWITCHES NOT INCLUDED / MUST PICK ONE FROM Pg. 9

5 - Sensor Output/ Power 112

This five switch adapter package is a perfect option for a stationary setting. The wall plug allows you to plug directly into power.

Part# SO-5PR+112 \$225.00



SWITCHES NOT INCLUDED / MUST PICK ONE FROM Pg. 9

5-Sensor Output/ Power 130


This is a portable system which will allow you to run the five switch adapter from a power wheelchair.

Part# SO-5PR+130 \$315.00
Chair Specific



Package Options 1, 2 & 3


These Packages Are Ready To Go!

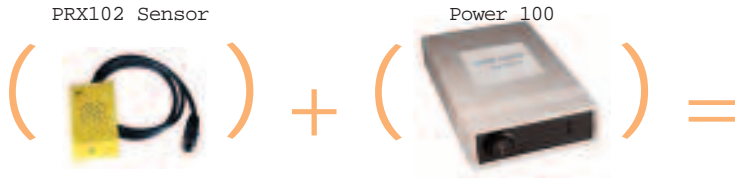

Package Option 1 - For Single Switch Scanners

Includes: Power 100, PRX103 proximity sensor used for single switch scanning option and Cable 1.


Part# PKG-01 \$400.00 (ACC device not included)

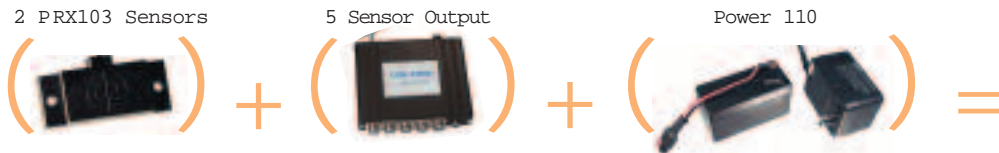

Package Option 2 - For Dual Switch Step Scanning

Includes: Power 110, two PRX103, Cable 4 with 5 sensor output. This is a great choice for dual switch scanning.

Part# PKG-02 \$680.00 (ACC device not included)

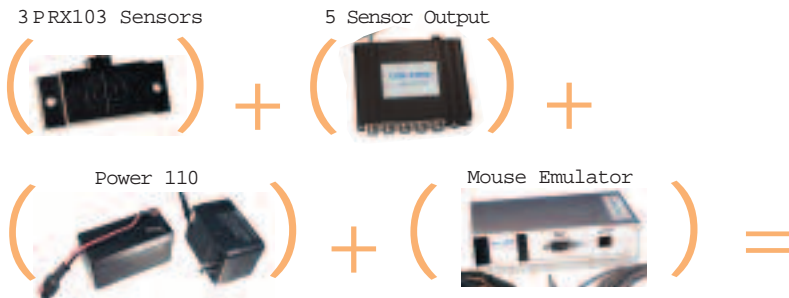

Option 3 - For Mouse Emulating

Includes: Power 110, three PRX103 proximity sensors, mouse emulator ME500 and Cable 4 with 5 sensor output. This is great access for Rolltalk and computer users.

Part# PKG-03 \$1,210.00 (ACC device not included)

comation Mechanical Switches



(light duty) MS-WL411

Wobble Switch

Wobble switches are activated by moving the rod in any direction. You can push it gently or whack it, the double spring helps to prevent damage to the switch. Select between a medium or light duty spring, which can be easily mounted for various positioning options.

Part# MS-WL411 \$65.00



MS-WM412
(medium duty)

Part# MS-WM412 \$65.00



MS-ML410

Micro Lite Switch

This light touch switch requires very minimal pressure and movement for activation. It can be easily mounted with velcro or permanently with a screw.

Part# MS-ML410 \$55.00



Egg Switch

The egg switches unique shape allows these switches to be mounted in close proximity to each other. The egg switch will activate with pressure on any part of the surface or on the edges. The cable comes out of the bottom and is routed through a trough so it is flush with the surface. There are threaded holes for screws to allow mounting or a ring option so the client can have the switch in their hand at all times. The egg switch comes in 5 primary colors.

Part# MS-400, MS-401, MS-402, MS-403, MS-404 \$37.00
(Black) (Red) (Yellow) (Blue) (Green)



Ring Mount

The ring mount option for the egg switch is a great mount for someone that needs to keep their switch with them. It attaches to the bottom of the egg and has a ring that slides on your finger. It comes in one adjustable fits all size to fit anyone.

Part# MS-RM405 \$17.00

(Egg switch not included)



Headrest Swing Away Bracket

The swing away arm attaches to any Stealth or Otto Bock style headrest. It is made ready to mount an egg switch or a sensor in a spot pad. Once you adjust it in place you can swing it away when you don't need it and just push it back into the exact same spot when you are ready.

Part# SM-HRSA \$240.00

(Headrest and egg switch not included)



The Button Switch

This amazing little switch is light pressure so it does not require much force to activate. It is shown here with three different mounting options, on a straw for use as a lip switch, on a T-handle and on a straight handle. It can be used as a reset option or mouse click with any 1/8 mono port.

Part# MS-BSS420 Straw Mount \$42.00 (Switch included)

Part# MS-BST422 T-Handle \$42.00 (Switch included)

Part# MS-BSI425 Straight Handle \$42.00 (Switch included)

comation Sip n Puff Switches




4 Function Sip "N" Puff Switch

Finally, a four-function switch in one. This revolutionary sip and puff switch and is great for use while in bed or while at a computer. In conjunction with a mouse emulator any hard puff would move the mouse up, soft puff right, hard sip down, soft sip left and with dwell or another switch site for click you are good to go.

(Requires power from power options page)

Part# MS-4FS454 \$360.00 * Requires breath tube Kit see below.

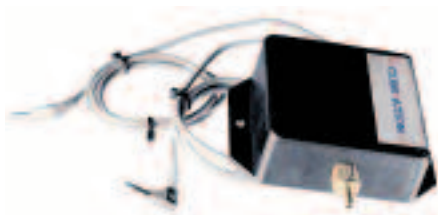

2 Function Sip "N" Puff Switch

This switch allows the user to have 2 switch activations by using sipping and puffing in a breath tube. These 2 switch activations are a great way for a client to left click and right click a mouse. There are 2 mono jack plugs so you can virtually plug it in to any device.

Part# MS-2FS450 \$180.00 * Requires breath tube Kit see below.

Breath Tube Kits

(All Breath Tube Kits come with a filter, oxygen tubing, and a straw.)



(back cane mount shown)

Back Cane Mount Breath Tube Kit Headrest Mount Breath Tube Kit

Both kits use flexible snaplock tubing which allows the breath tube to be positioned as necessary.

Part# SP-BCM451 \$237.00 *(back cane mount)*

Part# SP-HRM452 \$237.00 *(headrest mount)*



Ear Mount Breath Tube Kit

This headset is a perfect system for a client that needs to have the breath tube with them wherever they move. The headset mounts around the ears and the back of your head. Therefore, wherever you move your head the breath tube will be with you.

Part# SP-EM453 \$237.00



Bed Mount Breath Tube Kit

This flexible breath tube kit is designed for use while a person is in bed. Clamps to a bedside table or hospital bed rail.

Part# SP-BM455 \$237.00



Stealth Ultra

This headrest is the perfect combination of positioning and switch mounting. The sub occipital component will help a client to find midline. The two lateral swing away components are ideal for sensor mounting and you can also add a sensor in the backpad. It is a perfect match with the ME500 or ME503/5.

Part# SM-SU200 \$760.00



Stealth Elite

The elite is simple, beautiful and is great for someone that does not need a midline cue. The two lateral swing away components are perfect for switch mounting and if you add a sensor in the backpad it is a perfect match with the ME500 or ME503/5 for computer access.

Part# SM-SE200 \$590.00



Back Pack Mounting Bag

This bag will neatly hold all your electronics and all of the cables.

Part# SM-BG250 \$250.00



Table Clamp

This versatile mount clamps directly to a table or a desk edge. It has three joints so you can adjust it almost anywhere. This platform attachment allows you to attach any sensor with a power option to it.

Part# SM-TC220 \$155.00



Sensor Spot Pad Mounts

Spot pads are a great place to mount the PRX103. Embedded in this pad there are no sharp edges and the sensor range is still touch to 1/4 inch above the pad surface. They are easily mounted with velcro or by threaded holes on the back. (without sensors)

Part# SM-SP \$69.00 / Part# SM-MP \$69.00 / Part# SM-LP \$69.00



Swing Away Mount

The swing away mount attaches directly to a wheelchair armrest. With the pull of a lever the arm swings away to remove when not in use. The sensor and pad must be ordered separately.

Part# SM-SAM230 \$162.00



Headrest Swing Away Bracket

The swing away arm attaches to any Stealth or Otto Bock style headrest. It is made ready to mount an egg switch or a sensor in a spot pad. Once you adjust it in place you can swing it away when you don't need it and just push it back into the exact same spot when you are ready.

Part# SM-HRSA240 \$240.00 (Headrest not included)

comation Cable Connections CC, Computer Access & Mouse Emulators



Cable 2 Part# CC-602 \$65.00

This cable is a perfect match with a two switch scanning system. One end will have a D9 plug that can plug into the Power 5 or a power wheelchairs auxiliary output module. There are two mono jack plugs on the other end that will plug directly into the device. It is required that we know which device that you are plugging into because not all devices require the same input.



Cable 4 Part# CC-604 \$65.00

This cable gives you the most flexibility of all. One end has a D9 connector and there are four mono jack plugs on the other. It works great with joystick mouse emulation. In conjunction with a power wheelchairs auxiliary control module or a 5 sensor output you can hook up almost anything.



Cable 3 Part# CC-603 \$65.00

Paring this cable with a three switch mouse emulator is a perfect combination. One end is a D9 connector and there are three mono jack plugs on the on the other end. It is required that we know which device that you are plugging into because not all devices require the same input.



Cable 4 Part# CC-DX604 \$65.00

This cable will work with any Dynavox 3100 series devices. A standard cable will not function with the Dynavox 3100. If you want to use the driver control system as an access method, this cable is the one you need. The DX cable will always have a blue band so you can recognize it is for the Dynavox.

Lap Top Power Converter

This is a wonderful addition to a power wheelchair. A client will be able to use their laptop whenever they need it. No more dead battery. Simply plug the power converter into the power wheelchair and the computers cigarette adapter into the converter and the computer will stay charged all day. The converter does not drain power chair batteries so you just charge the chair in your normal routine.



Part# CA-809 \$525.00

*Requires wheelchair power source, Power 130. Specify wheelchair type. Specify computer amperage.

500 Series- 3 Or 5 Switch Mouse Emulator

The mouse emulator is a great way to access a computer. The ME500 can be set up as a three or five switch mouse. This version of the mouse emulator must be directly connected to the computer through the USB port and into the switch interface through a cable. This will enable the user to use three or five switches as a mouse on the computer. The three switch setting will toggle directionally. One switch will toggle up and down, the second switch will toggle left and right, the third switch will be click, double click, and drag. In the five switch version each switch is directional. One switch is up the second switch is down, the third is right, the fourth is left and the fifth will be click, double click, and drag. All of the mouse functions you need in one little box.



Part# ME-500 \$365.00

3 Switch Wireless Mouse Emulator

The MEIR 550 three switch wireless mouse emulator. One switch will toggle the mouse up and down the second will toggle the mouse right and left, the third switch will be click, double click, and drag. The receiver attaches through a USB port the computer can still have a standard mouse attached. The transmitter attaches to the auxiliary port of the power wheelchair. This wireless feature will enable the client to drive right up to the computer and start using it.

Part# ME-IR503 \$840.00

Switch Wireless Mouse Emulator

This five switch version is perfect for someone who uses a joystick to access their power chair or an external joystick attached to the computer. The receiver attaches through a USB port. The transmitter attaches to the auxiliary port of the power wheelchair. The client will have to have a dedicated switch for mouse click. The wireless feature on this mouse emulator will give the client the freedom to drive right up to the computer and have immediate access, without cables or wires! No more waiting to be plugged in.

Part# ME-IR505 \$840.00



On the cover

Salome drew this picture of "The Biloxi Lighthouse". It took her about three months to complete. Her mother stated, "We had to fold the paper into one inch squares. She never saw the finished product until she was completely through". Salome got the idea for her drawing from a postcard Mama found at a drug store. We think she did a great job don't you?



comation Fun N Games



Play Station 2

Play Pak

Now you can operate Sony's Playstation 2 with your drive control system from your power chair via the ECU unit on the chair. There are up to 12 mono ports that are available to control all the playing functions of the Playstation joystick on the interface. An added bonus is that the modified joystick functions can still be used if assistance is needed. Mechanical switches can be plugged into the interface as well. ECU unit and CC604 cable not included.

Part# FG-PSM700 \$150.00



Remote Control Toy

With this modification you can operate a remote control toy through your drive control system on your power chair. An ECU module or an auxiliary output module is required. You will need to send the toy in to be modified.

Not all toys can be modified, contact Comation first.

Part# FG-RFM710 \$150.00

WHEN DINOBOT ATTACKS

Mitchell and Devon
Thanks to the Rolltalk, Mitchell and Devon get to experience turn-taking and the thrill of cooperative play!
What better way to learn switch access than to play with "off the shelf" infrared toys like their peers. Mitchell and Devon definitely agree!!





comationTM

[changing Lives With Technology]

P.O. Box 255 125 Spur 191 Suite D Spicewood, Texas 78669
phone: 1(830)798.1914 fax: 1(830)693.6021 e-mail:comationinc@yahoo.com

COM ATION was founded to provide access products for expanding an individual s capabilities. Our product line reflects a growing number of components which will enhance interfacing computer, communication, mobility, and environmental technologies. Our focus is directed to SPECIFIC individual needs... not creating a standard method of access. Our knowledgeable and friendly customer service staff can discuss options of opportunity with you.